

First of all,

If you must have an Android version over 4.2.2


This App was created to have a direct view of the azure data without accessing the website. This is an app oriented to TABLETS not to phones however if your smart phone has a big screen this app could be used also in it.

Install



- Download the apk (<https://github.com/arbox0/NightClock/releases>) or download the source code and compile your own apk.
- Use your file browser on your phone to install the apk.

Config


The first time that you execute NightClock you will see this screen:

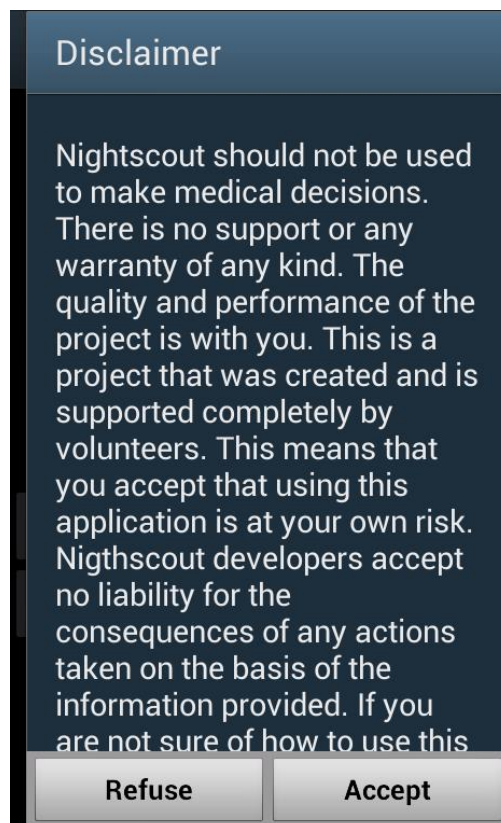


Please read carefully and if you agree, press Accept button, otherwise press Refuse button (this last option will kill the application).

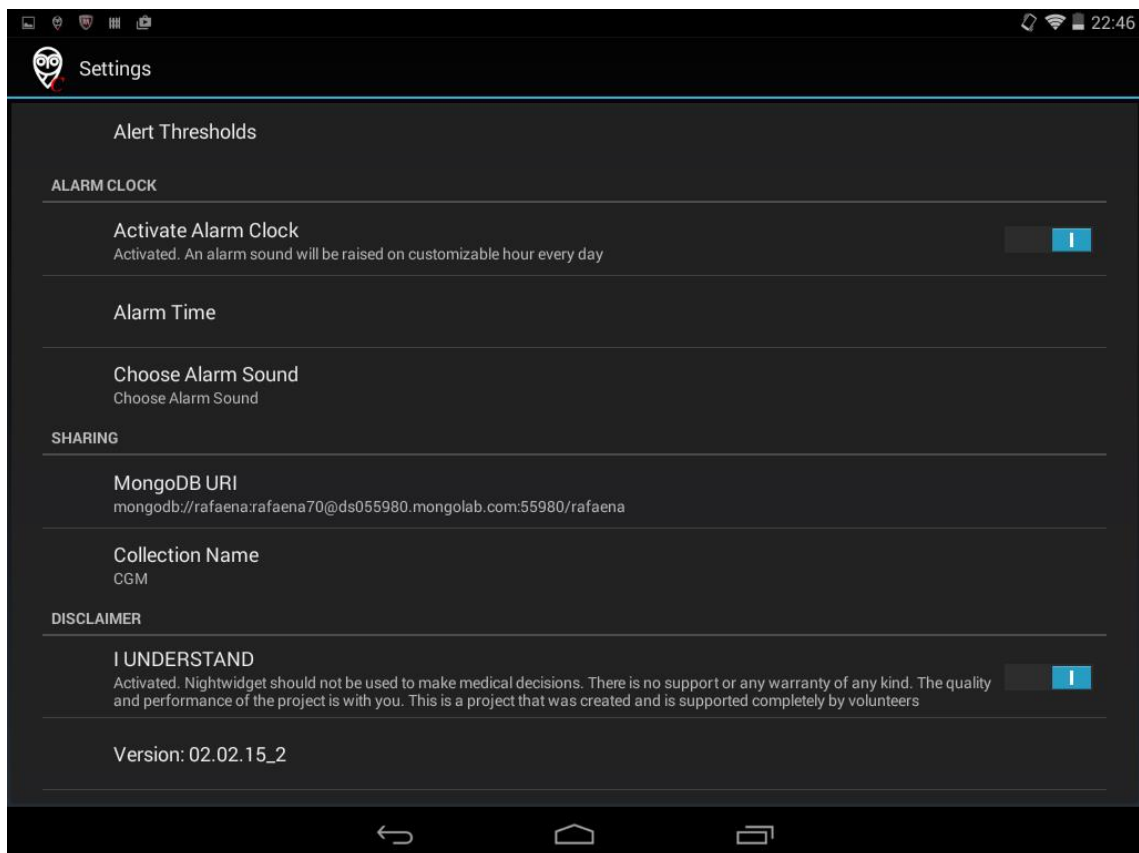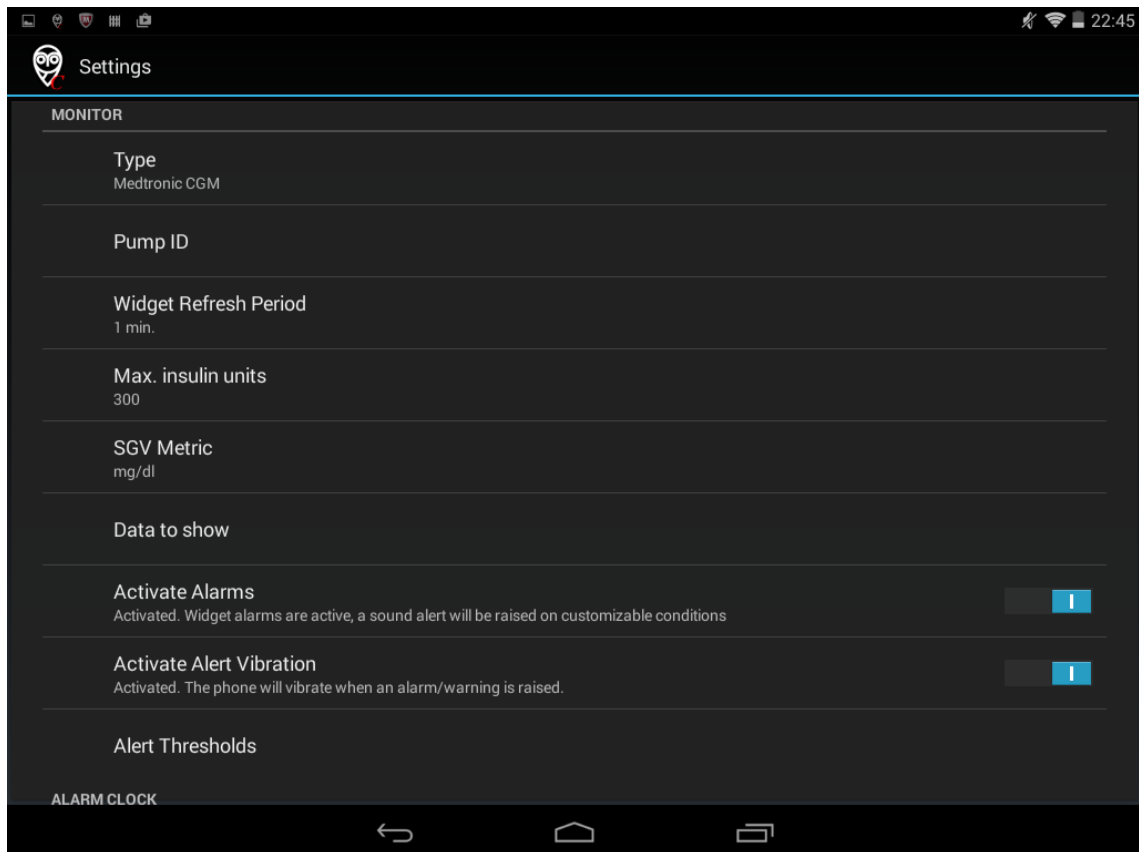

1. **Type.** Select your device type:
 - a. Medtronic
 - b. Dexcom
2. **Medtronic CGM ID.** ONLY if you have selected Medtronic Type. You must set the ID of your Medtronic pump to get access to Medtronic additional info.
 - a. **Widget Refresh Period.** The nightclock will be refreshed each: 1 Min., 2 Min., 3 Min., 4 Min., 5 Min., 10 Min., 15 Min., 20 Min., 30 Min., 60 Min.
3. **Max Insulin Units.** ONLY if you have selected Medtronic Type. You must set capacity of insulin units in the pump.
4. **Data to Show.** This preference shows a menu to select which options will be shown on your screen if data is available on your MongoDB.

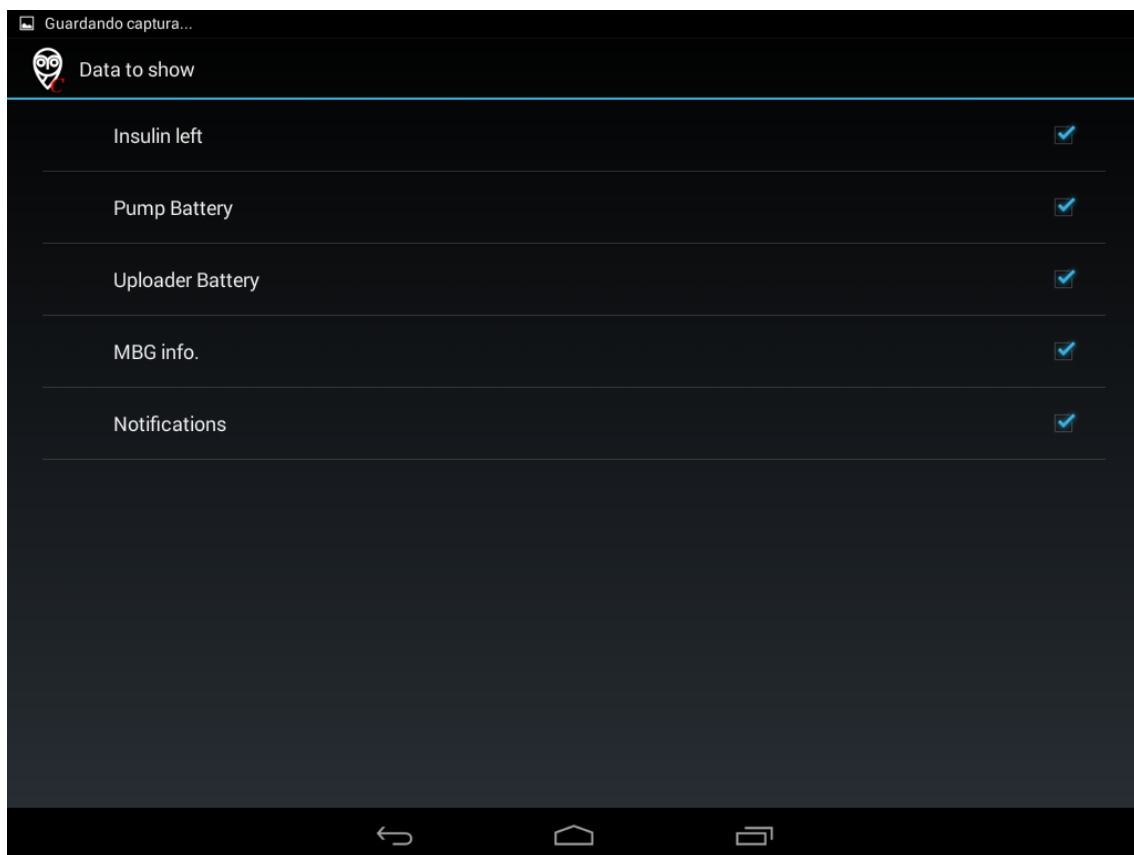


5. **Activate Alarms.**
 - a. **Off.** No alarms will be raised
 - b. **On.** Clock alarms are active, a sound alert will be raised on customizable conditions.
6. **Alert Thresholds (Submenu).** Enabled if Activate Alarms is On, Disabled otherwise.

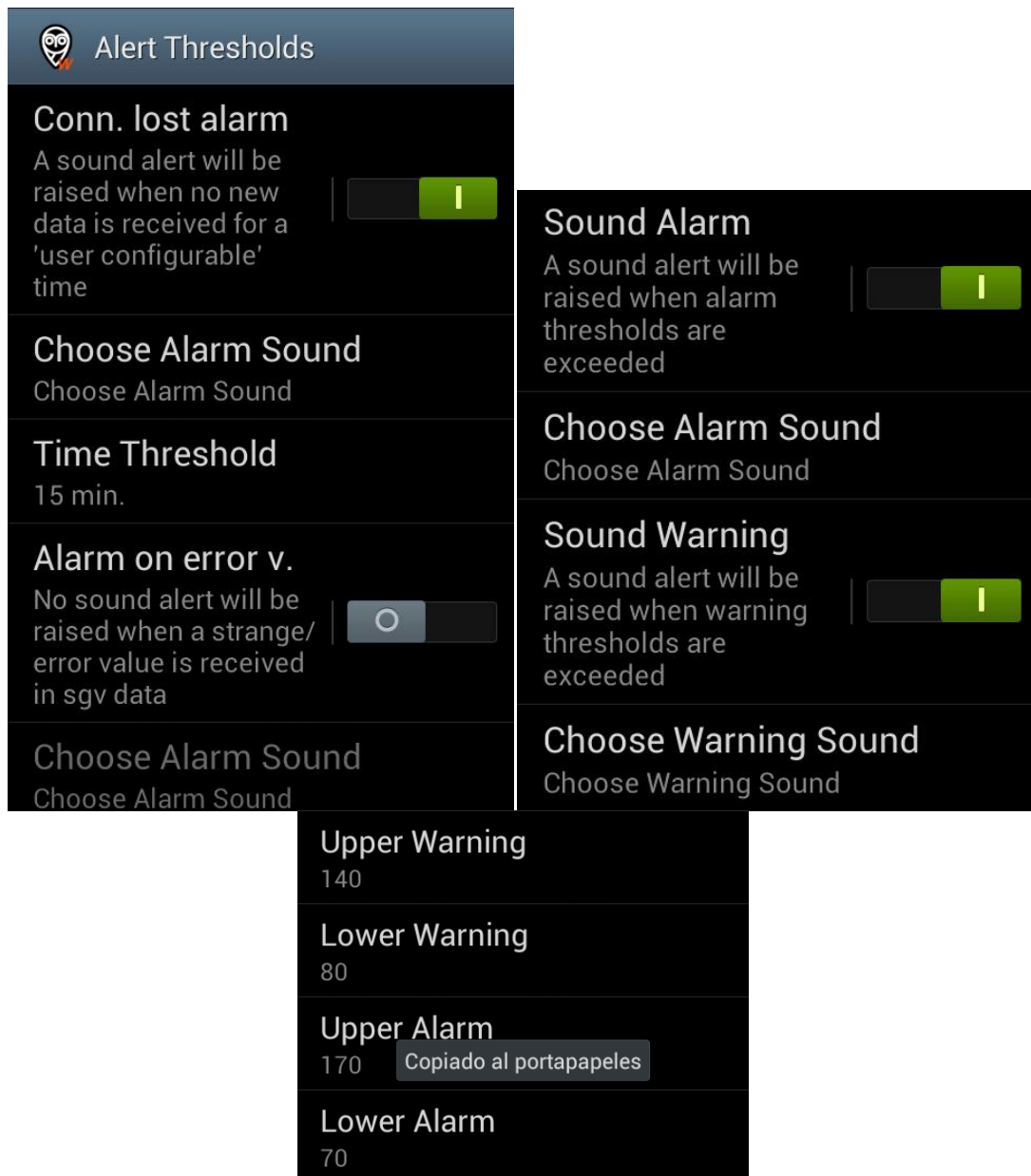


- a. **Connection Lost Alarm.**
 - i. **On.** A sound alert will be raised when no new date is received for a "user configurable" time.
 - ii. **Off.** No alert will be raised.
- b. **Choose Alarm Sound.** Set the alarm sound track to be played as alarm.
- c. **Time Threshold.** Time after which the nightclock will raise the connection lost alarm.
- d. **Alarm on Error value.**
 - i. **On.** A sound alert will be raised when a strange/error value is received as sgV data.
 - ii. **Off.** No alert will be raised.
- e. **Choose Alarm Sound.** Set the alarm sound track to be played as alarm.
- f. **Sound Alarm.**
 - i. **On.** A sound alert will be raised when alarm thresholds are exceeded.

- ii. **Off.** No alert will be raised.
 - g. **Choose Alarm Sound.** Set the alarm sound track to be played as alarm.
 - h. **Sound Warning.**
 - i. **On.** A sound alert will be raised when warning thresholds are exceeded.
 - ii. **Off.** No alert will be raised.
 - i. **Choose Warning Sound.** Set the alarm sound track to be played as warning.
 - j. **Upper Warning.** Upper Warning Threshold.
 - k. **Lower Warning.** Upper Warning Threshold.
 - l. **Upper Alarm.** Upper Alarm Threshold.
 - m. **Lower Alarm.** Lower Alarm Threshold.
7. **MongoDB URI.** Set the address to your Mongo database (It must be the same address used in Nightscout uploader).
 8. **Collection Name.** Set the sgv entries collection name in the mongo data base. (It must be the same collection name used in Nightscout uploader)
 9. **Disclaimer (I Understand).**
 - a. **On.** You must accept to activate the nightclock functionalities.
 - b. **Off.** The nightclock is turned off.

NightClock View

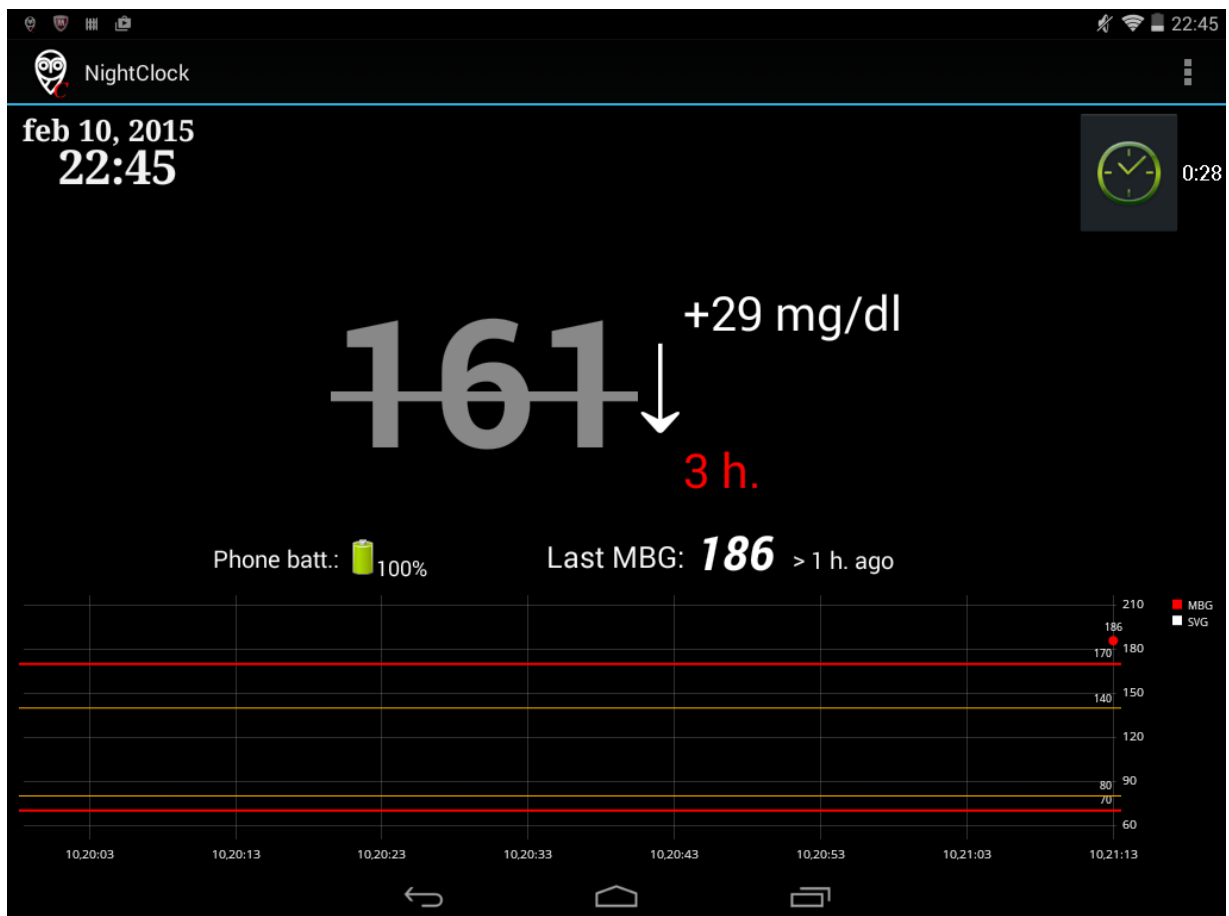

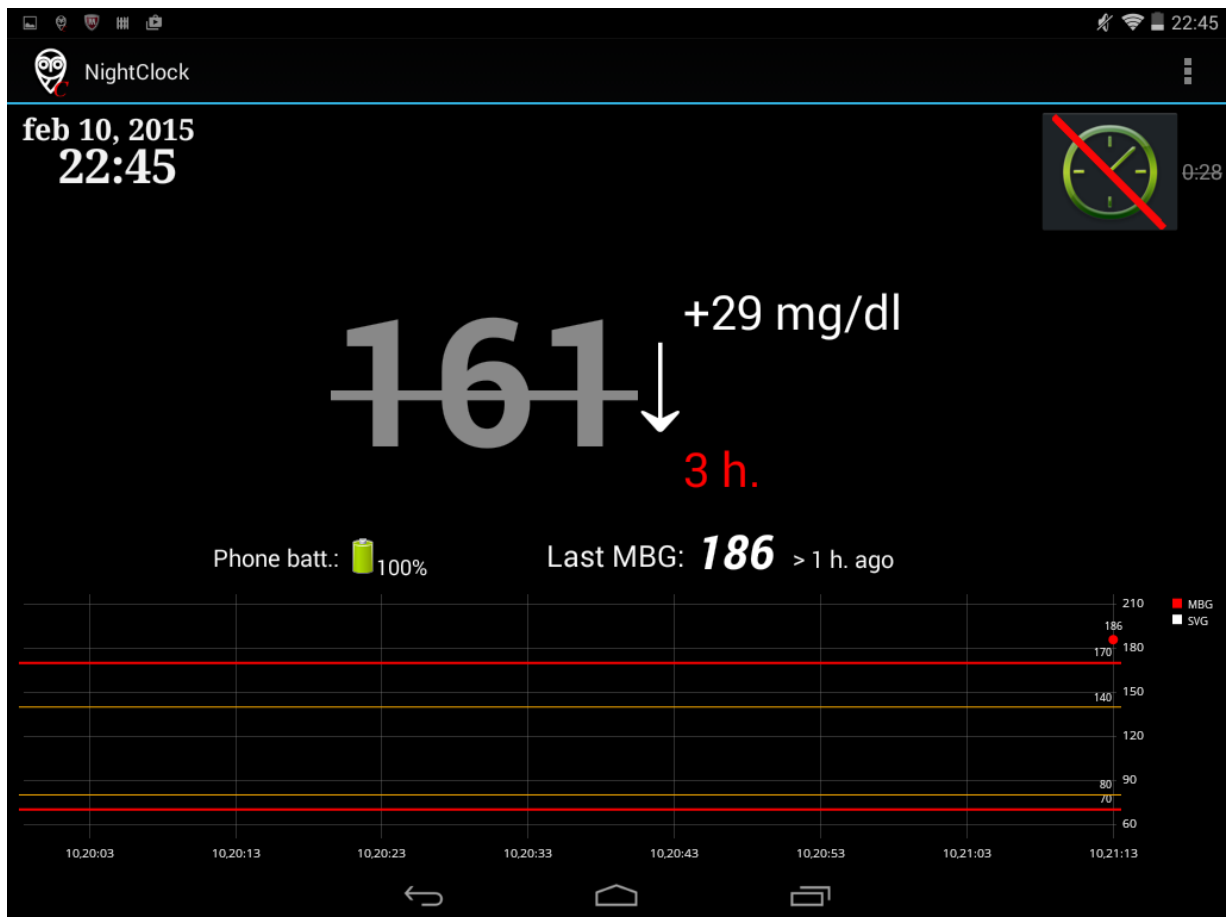

If you see:

- 125?! → Last calibration has failed, value received 125
- 125! → not ideal calibration but good enough, value received 125
- 125? → last calibration is > 12 hours old, value received 125
- 125* → is calibrating, value received 125
- 125+ → is calibrating second stage, value received 125
- NC → not calibrated
- DB? → database error??